




SOUTHERNIMPLANTS®

Innovative Treatment Solutions

BHORA GUIDED SURGERY KIT


Precision Guided. Universally Compatible.

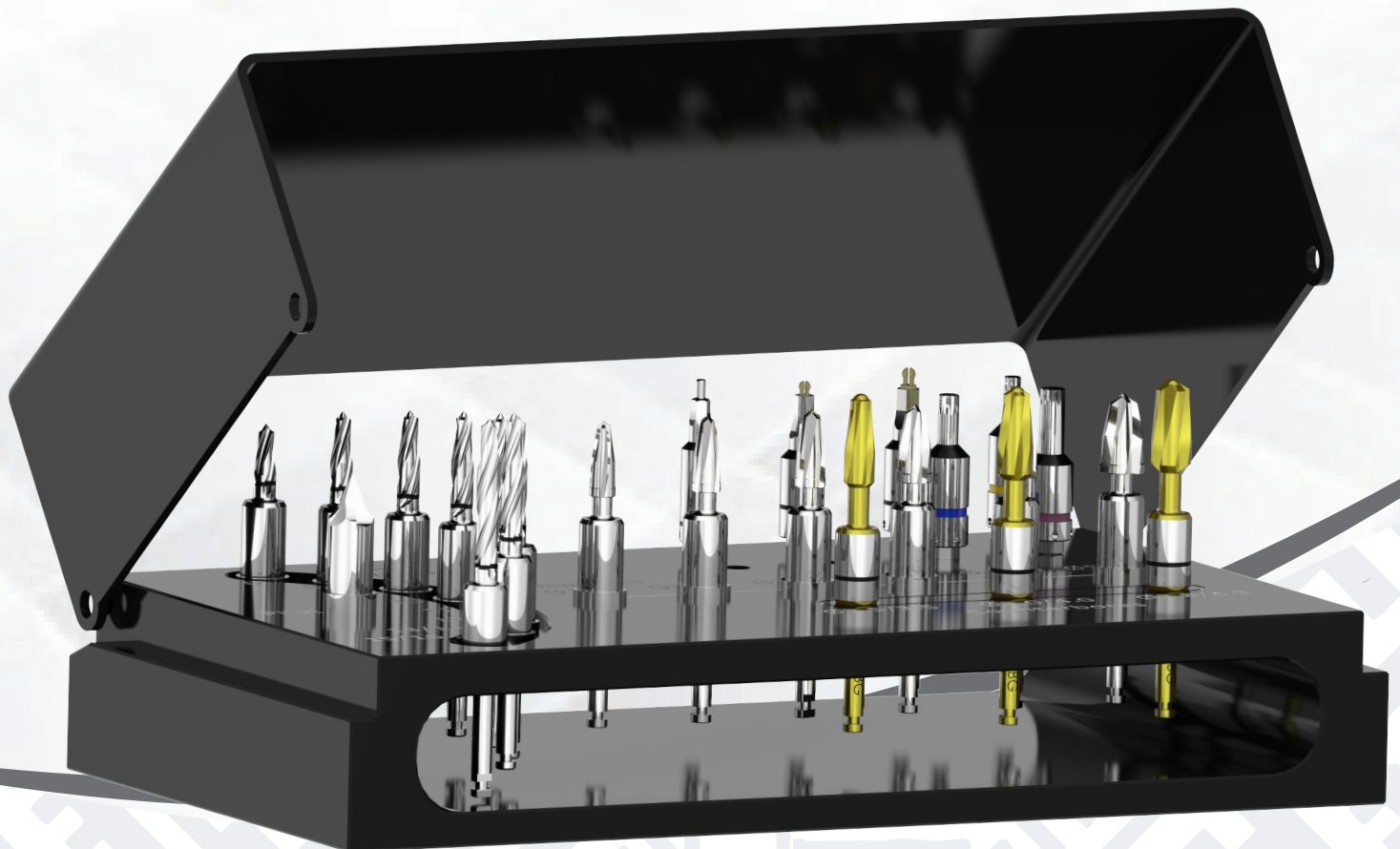


BHORA Guided Surgery Kit

Designed for clinicians.
Built for predictable
simplified guided workflow.

Meet the
**SIMPLIFIED UNIVERSAL
GUIDED SOLUTION**
for placing Southern
tapered implants.
($\varnothing$ 3-5 mm, 8-16 mm)

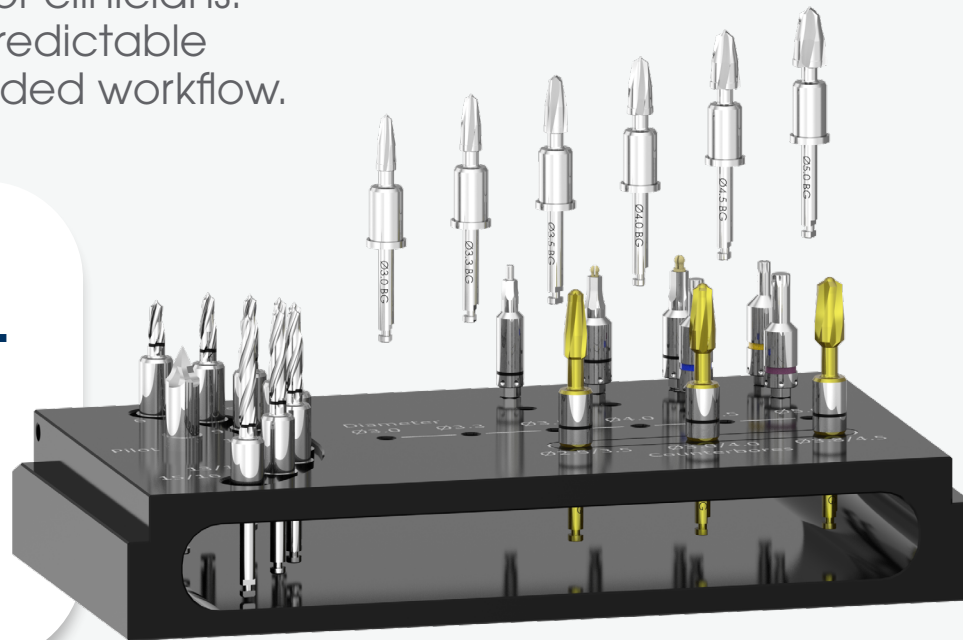

One offset. Total control.
**10.5 mm sleeve offset with
built-in stops.**

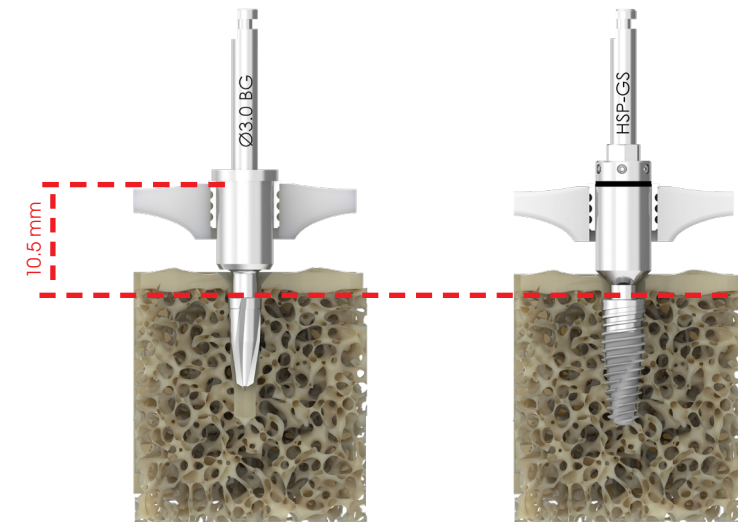

Note: Lip on the guide sleeve is 0.25 mm thick. Offset should be taken from the top of the sleeve.

4 step
guided workflow.

Step 1

Pilot Drill



Initiate the osteotomy
by perforating
and profiling the
cortical bone.

Step 2


Depth Drills



Incrementally increase
the osteotomy
depth to the
planned final length.

Step 3


Diameter Drills

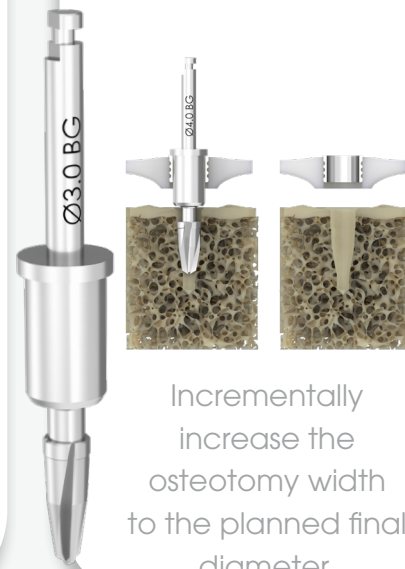

Incrementally
increase the
osteotomy width
to the planned final
diameter.

Step 4

**Wide Apex
Counterbore**
Optional

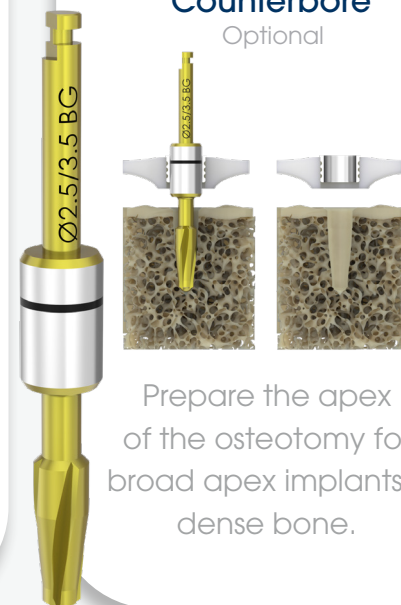


Prepare the apex
of the osteotomy for
broad apex implants in
dense bone.

Fewer drills.
Smarter design.


- ✓ Reduced inventory.
- ✓ More precise surgeries.
- ✓ A more efficient guided workflow.

Drill Protocol Calculator
Scan QR code to access our
drill protocol calculator.

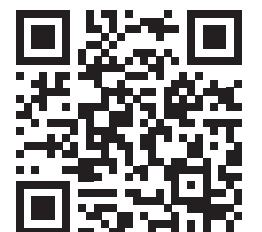

Order Guide

I-BG-MT

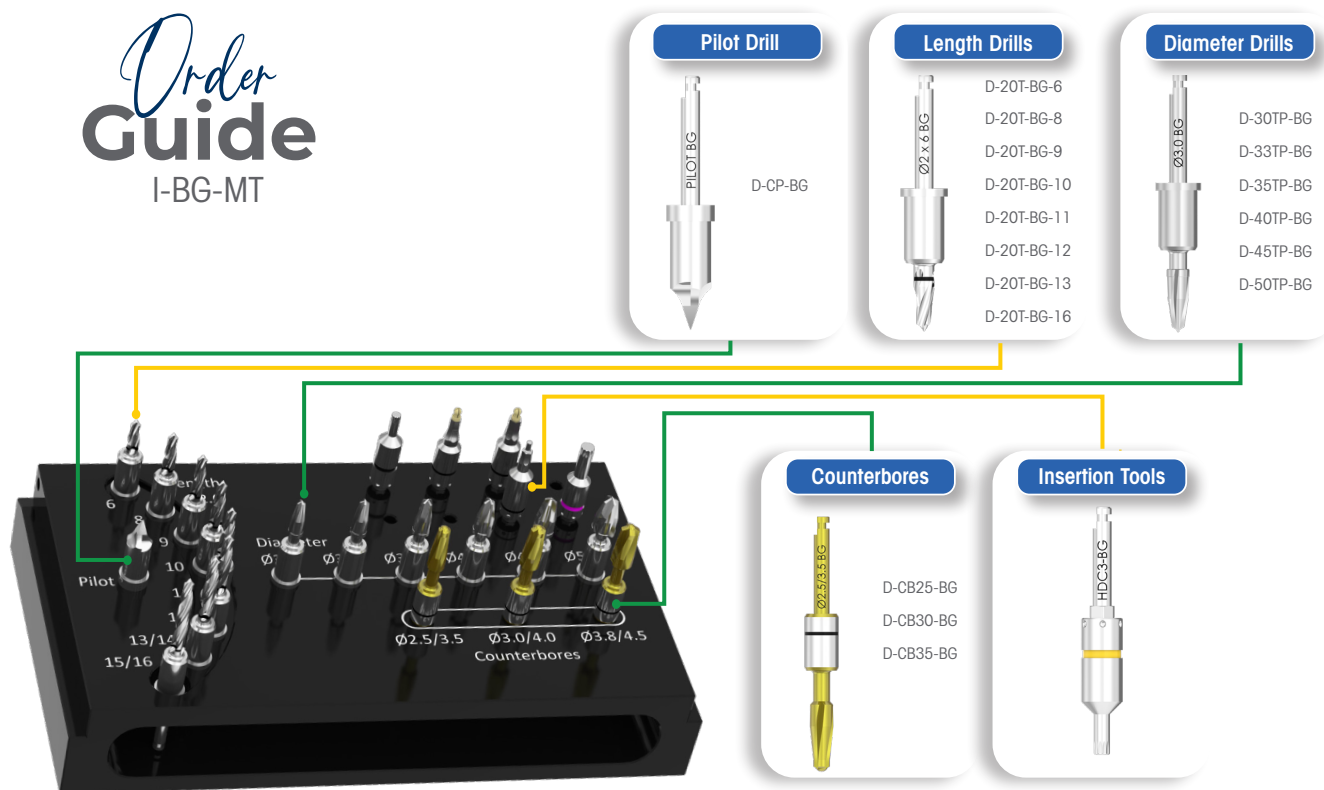

Standard Guided Sleeves Ø6.2 (outer) / Ø5.2 (inner)

Closed Sleeves



I-GS-6252-L4



I-GSF-6252-L4
anti-rotational flat
I-GSF-6252-L4-P5
5 pack



I-GSF-6252-L6
anti-rotational flat



I-GS-6252-L9



I-GSC-6252-L4



I-GSCF-6252-L6
anti-rotational flat



I-GSC-6252-L9

Guide Fixation Components



I-GS-6252-L4

Drill Ø1.2 x 15 mm



I-D12T-GP

Guide Fixation Pin



I-GS-2513-L6

Guide Fixation Pin Sleeve

Fully Guided Implant Placement Tools

DC (Deep Conical)

I-HDC3-BG	Ø3.0 DC	Straight
I-HDC4-BG	Ø4.0 DC	Straight
I-HDC5-BG	Ø5.0 DC	Straight

SP1 (Single Platform)

I-HSP-GS	Straight
----------	----------

Internal Hex (M-Series) and PROVATA®

I-H3M-BG	Ø3.0 PROVATA®	Straight
I-HM-BG		Straight
I-H-PRO12D-GS		Co-Axis®
Internal Octagon (IT)		
I-HITS-GS		Straight

TRI-NEX®

I-HLH-35GS	Ø3.5 TRI-NEX®	Straight
I-HLH-43GS	Ø4.3 TRI-NEX®	Straight
I-HLH-50GS	Ø5.0 TRI-NEX®	Straight
I-L43-12D-GS	Ø4.3 TRI-NEX®	Co-Axis®
I-L50-12D-GS	Ø5.0 TRI-NEX®	Co-Axis®



South Africa – Headquarters 1 Albert Road, Irene, RSA
T: +27-12-667-1046 | E: info@southernimplants.com

EC REP

Southern Implants Europe AB Holmgatan 30, S-791 71, Falun, Sweden
E: ecrep@southernimplants.com

Subsidiaries

Australia
Southern Implants Australia
T: +61-8-9466-2627
E: info@southernimplants.com.au

Spain and Portugal
Southern Implants Ibérica
T: +34-935-053-507
E: info@southernimplants.es

United Kingdom and Ireland
Southern Implants UK
T: +44-20-805-94490
E: info@southernimplants.co.uk

USA and Canada
Southern Implants North America Inc.
T: +1-561-472-0990
E: customercare@southernimplants.com

Southern Implants® are distributed world-wide, please visit southernimplants.com for a list of Distributors.
For the latest and full assortment of Southern Implants' products, visit our current resources library at southernimplants.com.