ScanLadder

Placement Instructions



1 Choose Scan Ladders & Insert Screws

Assess the space between each multi unit abutment in the dental arch. Insert the screw and hold with the screwdriver ready for insertion.




2 Place and Screw Scan Ladder to MUA

Next, insert the Scan Ladder and with one hand hold the scan body into position whilst the other screws the scan body into place using the screw driver included in the kit. Use the screw torque tightening force recommended by the manufacturer.



3 Centrally position the Scan Ladders

Each Scan Ladder should be positioned with the aim to centrally position the shaped end of the scan body OR follow the arch as is possible in the patient case. Choose the length accordingly and vrify seating on a radiograph.

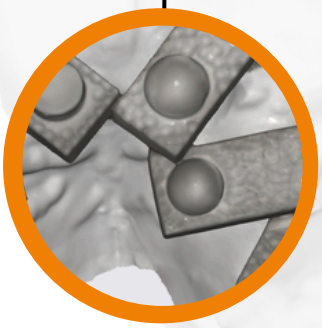

4 Avoid interferences & compression

Once youve placed each scan ladder, make sure that each scan body does not compress the adjacent scan body through rotation or overlapping. Overlapping is ok as long as it doesnt impede seating.
Check passivity with floss if needed.

Placement Examples

