
ScanLadder

Tips & Recommendations



1 The scanner tip shouldnt touch

Make sure to avoid bringing the scanner too close to the scan bodys or indeed knocking the scan body out of place. If you do knock a scan ladder loose, retighten and cut out that from the scan. Restart the scan from the central position for that scan body.

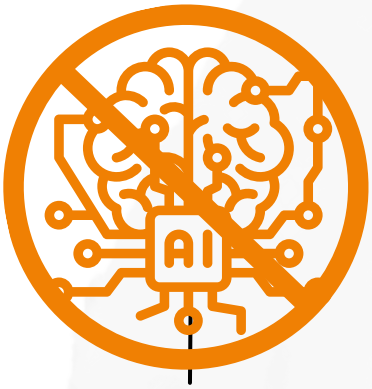


2 Keep scanning smoothly & Turn Off AI

As you scan following the path described in the “Scanning Path Instructions” document, make sure that the scan motion is smooth and not jerky. If your scanner has any AI functions, turn these off.



3 If you lose track, go back to the centre

If you lose track, do not go back to a cylinder, go back to the centralised position of the end shapes or to the arch ridge at a 45 degree angle with a characteristic and random surface.

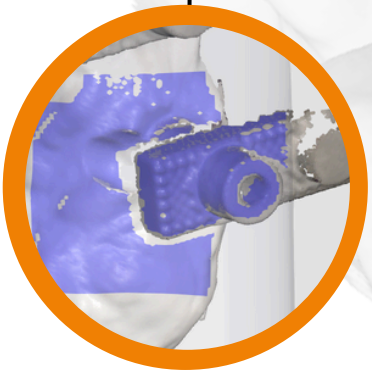

4 Dry the field and avoid gaps

Lower arches and bloody fields are conventionally more difficult. Dry the field and the scan body random surface, clearing any blood and saliva as much as possible with good aspiration control. Make sure that scan bodys are adjacent or focused near to each other to avoid gaps in a continuous scan.